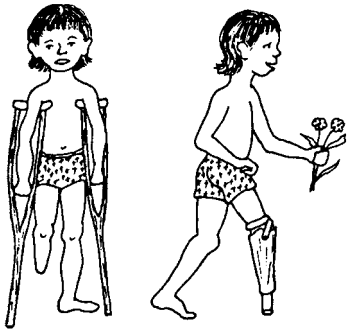


Artificial Legs

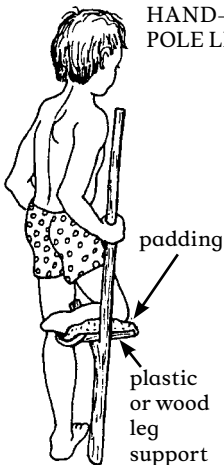
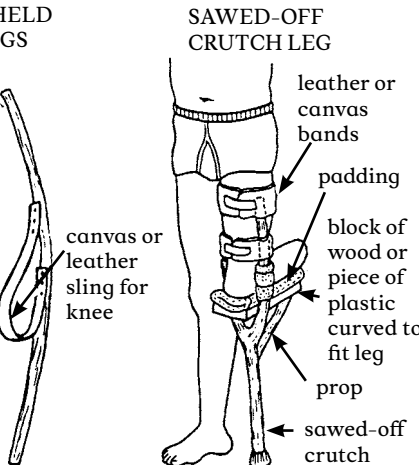

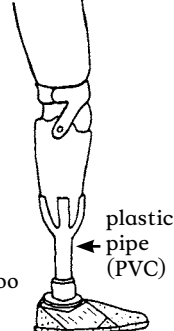
Artificial legs can be (and often are) made at home or in local shops. How well they work and how natural they look depend on many things, including costs, skills, and materials available.



BELOW THE KNEE

The most common leg amputation is below the knee. A leg that has been amputated halfway between the knee and ankle works best for walking with an artificial limb. Here are some examples of artificial limbs, from simple to more complex.



Even a simple artificial limb can make a big difference.


WITH THE KNEE BENT		WITH THE KNEE STRAIGHT	
HAND-HELD POLE LEGS	SAWED-OFF CRUTCH LEG	BAMBOO AND PLASTER LEG (see p. 628)	PLASTIC PIPE LEG WITH FOOT (see p. 632)
			
<p>CAUTION: Limbs like these 3 are quick and easy to make, but they cause knee contractures. As a result, the knee cannot be easily straightened to fit a better, more useful limb. Bent-knee limbs should only be for temporary or emergency use. Do exercises every day to keep the knee straight and strong (see pp. 229 and 230).</p>		<p>These limbs are better because the knee has full range of motion. Walking is easier and more natural. However, the person's weight must be supported evenly over the entire stump, not only at the end of the stump (see p. 631).</p>	

Positions for FITTING A LIMB	BAD	DIFFICULT	GOOD
	will only work with a bent-knee limb		knee straightens completely

Exercises to strengthen and straighten the leg

From the time a leg has been amputated until a limb is fitted, daily exercises are needed to keep the hip and knee muscles strong and to avoid contractures. If weakness and contractures already exist, these should be corrected as much as possible before a limb is fitted. Exercises are discussed on pp. 229 and 230.

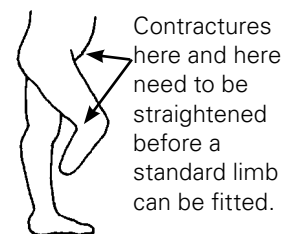AVOID CONTRACTURES



How soon can an artificial limb be fitted?

Children born without a foot or part of a leg (or legs) can be fitted with an artificial limb as early as 10 or 12 months of age.

A child whose foot has been cut off can and should be fitted with a temporary limb as soon as the wound has healed. However, be very careful not to injure or put any pressure on the new scars or end of the stump.



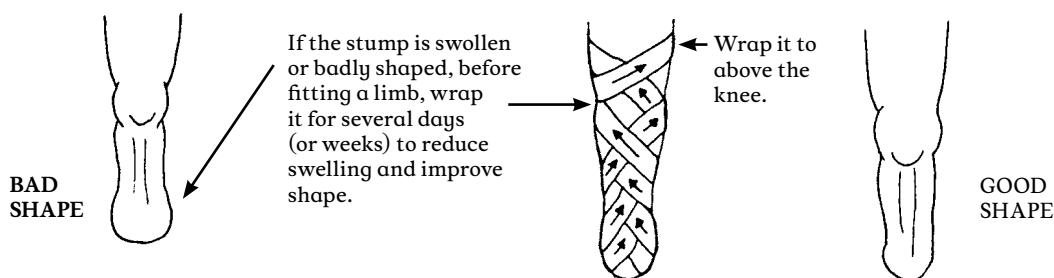
Note: On some children it may be difficult to fasten the limb firmly to the knee (the bones may not stick out enough). Straps to a waistband and even over the shoulder may be needed.

Temporary limbs—when to use them and why

Because a stump usually shrinks and changes shape in the first weeks after a limb is fitted, it is often wise first to fit a low-cost, temporary limb. This is especially true if the amputation is new or the stump is swollen. A better-looking, more permanent limb can be made after 4 to 6 weeks, or when swelling is gone.

Preparing the stump

In the first weeks or months after an amputation, the stump tends to swell up. The swelling may in time lead to a stump which is difficult to fit with an artificial limb. For this reason, it is important to wrap the stump with elastic bandage from the time the leg is cut off until a limb is fitted, or at least until there is no more sign of swelling. Instructions for wrapping the stump are on p. 228.




NOTE: When the person is not wearing the artificial leg, he should also wear an elastic bandage to control the stump shape.

The art of limb making

Making artificial limbs that fit and work well is both a science and an art. If possible, try to learn from a skilled limb maker. On-the-job training for even a few days can make a big difference.

Before starting to make an artificial limb, STUDY THE PERSON'S LEG.

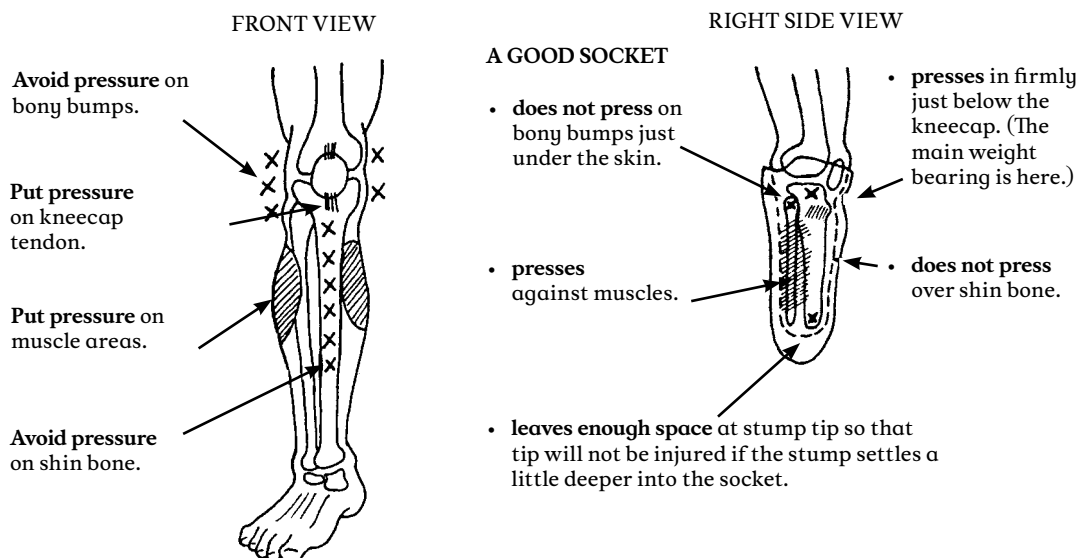
A good fit of the socket on the stump and at the knee is one of the most important—and difficult—parts of limb making. It helps to have an understanding of the bones and muscles in the leg



Before beginning, study the person's knee and stump carefully. Note the positions of the kneecap, the bony bumps on the sides of the knee, and the shin bone.

MEASUREMENTS YOU MAY NEED FOR A BELOW-KNEE LIMB

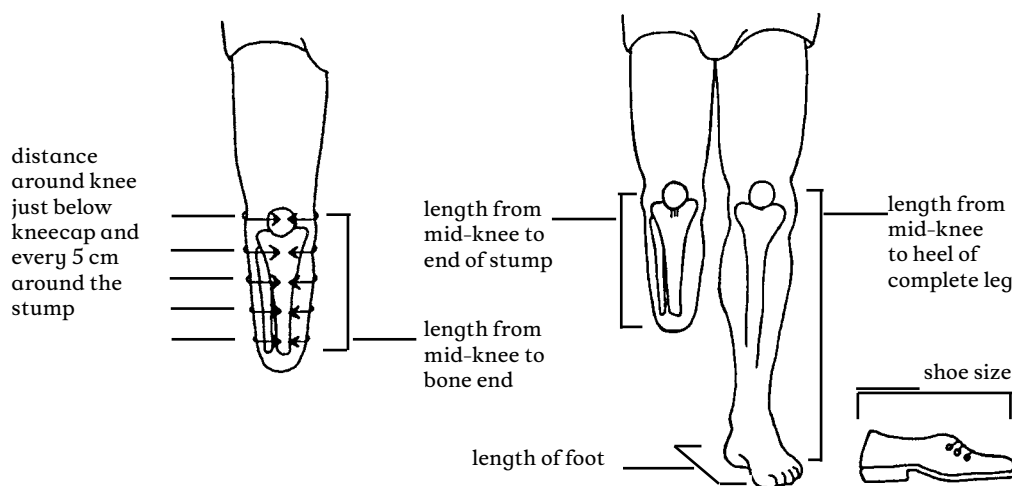(Copy this chart and use it to record your measurements.)



Note: For the plaster and bamboo limb, only lengthwise measurements are needed.

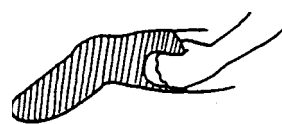
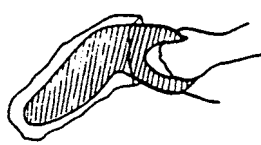
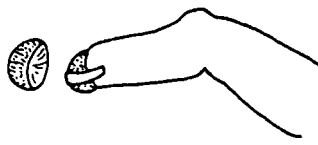
Note: The artificial limb should be the same length or just a little shorter than the other leg.

PLASTER AND BAMBOO BELOW-KNEE ARTIFICIAL LEG


This simple, low-cost leg was developed for refugees with amputations in Thailand by Operation Handicap International. It is most useful as a temporary limb for learning to walk. However, if the inner (plaster) part of the socket is made with waterproof glue, or is protected from getting wet, the leg can last for a long time.

Steps for making the plaster socket


1. Make a thick “cup” or “cap” of sponge or folded cloth and tape it over the end of the stump (to give it a little extra length).
2. Put a thick, tight-fitting stocking—without seams if possible—over the stump and knee. (Several thin layers of stocking can be used instead of a thick one.)
3. Put a thin plastic bag over the stocking.
4. Put a thin cotton stocking or stockinette over the plastic bag and also pull tight to avoid wrinkles.
5. With the stocking stretched tight, mark the important places with a grease pencil. The pencil marks will print onto the inside of the plaster cast when it is removed.



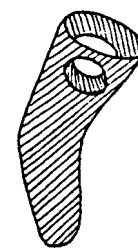A village rehabilitation worker fits a young man with a bamboo limb. (PROJIMO)




Putting holes near the top of the stocking makes it easy for the child to pull it up tight against the skin while the cast is put on.




Note: If you do not have a stocking, you can mold the socket directly on the stump. Shave any hair off the stump and cover it with vegetable oil (for example, coconut or palm oil).

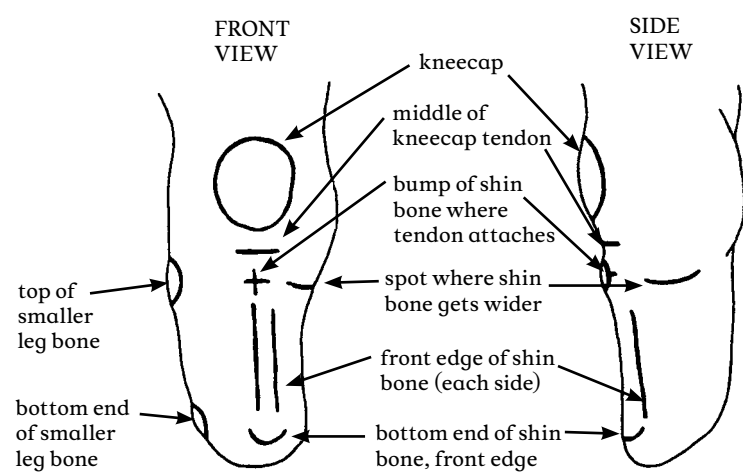Mark all these places:



6. Wrap the stump and knee with plaster bandage. Be very careful to put the bandage on evenly and smoothly. (Elastic plaster bandage works best, but is very costly. To reduce costs you can make your own plaster bandages for casting. See p. 569.)

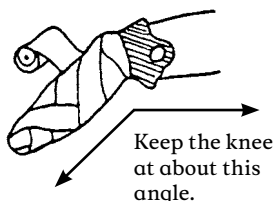Cast the stump with the knee slightly bent.



MAKING THE PLASTER WATER-RESISTANT

The plaster cast of the stump will become the inner layer of the socket of the bamboo limb. So it should be strong and waterproof. To make the cast stronger and water-resistant, wet the plaster bandage with glue instead of water.



Use a water-base glue that is water-resistant when it dries.

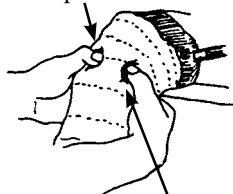
NOTE: If the plaster cast is to be used only as a mold for making a leather or resin socket, use water, not glue.

7. As the plaster dries, hold the stump firmly below the knee.



Press your thumbs on each side of the tendon just below the kneecap.

FRONT VIEW



With the heel of the hand, press gently on the muscles on both sides of the leg.

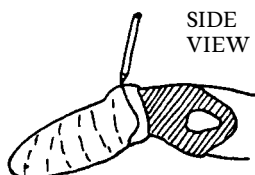
BACK VIEW



With the 2 middle fingers of each hand, press into the hollow behind the knee.

Hold the stump like this until the plaster is hard enough to hold its shape.

8. When the cast becomes hard, mark where to cut the top edge (see below).



SIDE VIEW

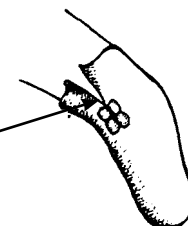


BACK VIEW

9. To remove the cast, roll the stocking over it. Put your hands over the pressure points (as shown above). Have the child wiggle the stump as you gently pull off the cast.




It may be necessary to cut the cast behind the knee, like this, to get it off.



10. Cut the cast along the line you drew.

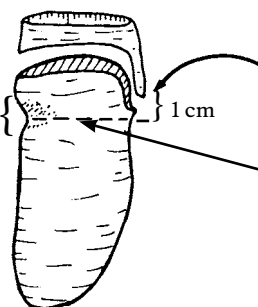
3 to 4 cm




The dip behind the knee should reach at most 1/2 cm above a line passing through the mid-point of the kneecap tendon.

Preparing the bamboo post

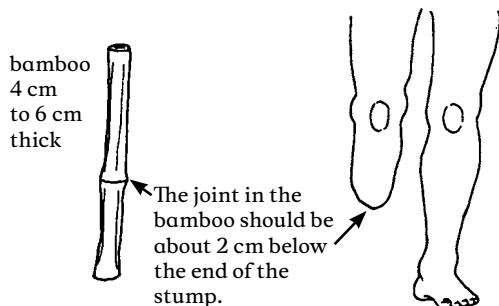
1. Select a piece of strong, green bamboo a little longer than the good leg from the knee to heel.



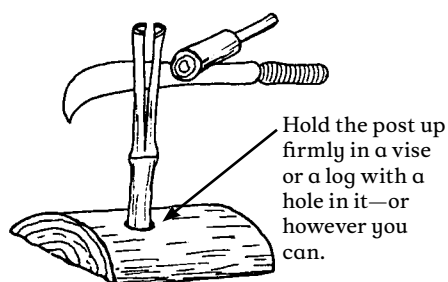
2. Split the bamboo to a little below the level of the stump end. Split into thin strips—each about 3/4 cm wide.



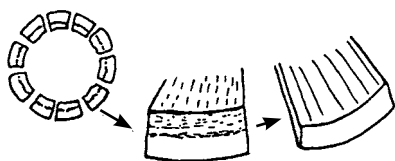
3. Remove the softer inner layer from each of the thin strips.



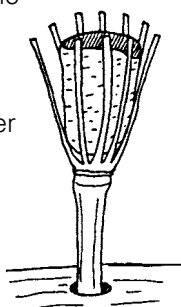
4. Spread the bamboo strips around the plaster socket.



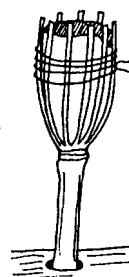5. Position the socket as correctly as you can. Then, with a thin wire, wrap the bamboo tightly against the socket.



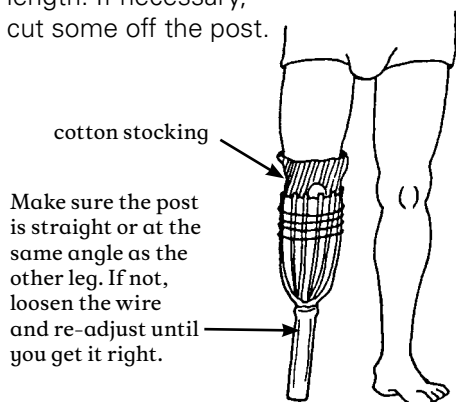
6. Put the limb on the stump and have the child stand on it. Check the length. If necessary, cut some off the post.




NOTE: If child uses a shoe or sandal, be sure to have her wear it when measuring the height of the limb.

7. After trimming the tops of the bamboo strips, cover the outside of the socket with several layers of glue, sawdust, and gauze bandage:

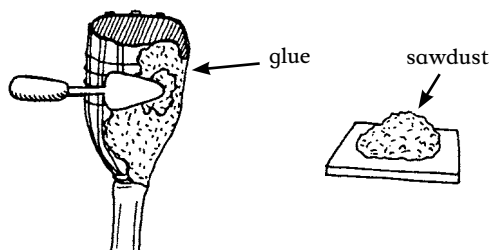- Brush on one layer of glue.
- Press sawdust on the glue (with gloves).
- Wrap tightly with gauze bandage. Let it dry.
- Repeat 5 or 6 times.




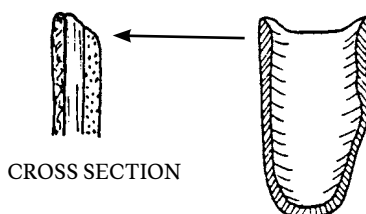8. Smooth the outside of the socket with sandpaper. Also smooth and round the inner edges at the top.




9. Make a cuff to hold the limb on. (If attached correctly, it should also help keep the knee from over-straightening.)

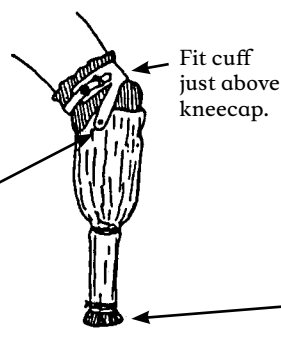Rivet cuff to socket slightly behind mid-line. (Mark rivet points with tabs pulled tight when knee is a little bent.)



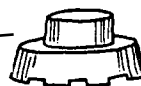
10. Make a rubber heel—a piece of thick truck tire works well. If you can, cut the tire so that a plug fits inside the bamboo. Cut off bamboo as much as the heel is thick. (Be sure to allow for height of sandal or shoe on other foot.)



Design for the knee cuff

Use strong leather and line it on inside with soft smooth leather. It should be a few cm longer than the distance around the knee.

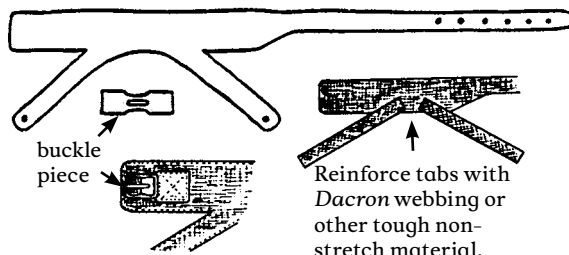Glue layers of leather and webbing together with rubber cement and sew at edges.



For other ways of holding the limb to the stump, see p. 626.

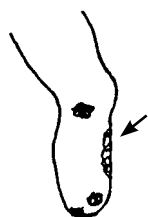
11. Have the child stand and walk on the limb for several minutes. Then remove it and look for sore spots on the child's skin, or signs of too much pressure. Check especially over bony spots. An area that looks pale when the limb is removed and then turns red or dark, is a sign of too much pressure.

12. Scrape shallow pits into the socket from the inside, at the points where it presses over bones. You may also need to build up around the area where pressure occurs.




To help you find the right points, it helps to mark bony spots with chalk or a grease pencil before casting. They will then be printed inside the socket (see p. 628) to show you the points that need scraping.


If stump presses on bottom of socket, you may need to build up. To build up these spots, dig a few pits into the socket surface so new material will grip better. Fill with a paste of fresh plaster mixed with glue.



bony spots

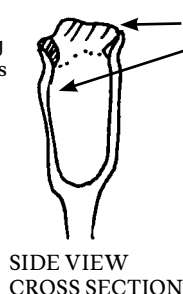

SIDE VIEW
CROSS SECTION

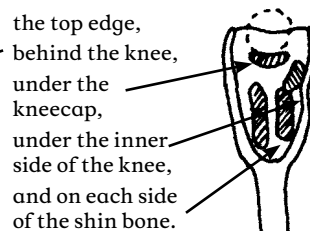

BACK VIEW
CROSS SECTION

If socket presses on the end of the bone, scrape out a hollow here.



During the first few weeks of using an artificial limb, the stump becomes smaller, and several changes in limb size may be needed. To save time, use a shorter bamboo post so that the plaster socket can be replaced several times with new, smaller ones.

Put the glue and sawdust only up to here.

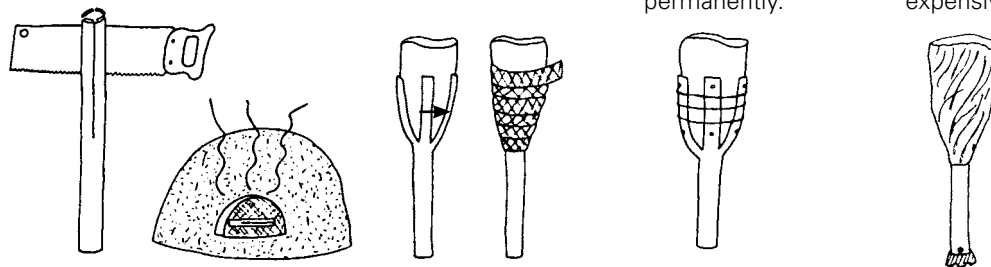


Hold the bamboo strips to the cast with wire. Cut the wire to change the plaster socket.

Artificial leg using PVC plastic pipe instead of bamboo

Where plastic PVC water pipe is available, it can be used instead of bamboo. Use a 3 cm (1 ¼ inch) thick-walled PVC.


1. Measure the tube the same as bamboo, and cut it to form 4 strips.
2. Heat the PVC in an oven until it gets a little soft.
3. Fit the hot PVC around the socket piece and wrap it tightly with a long strip of cloth or rubber until it cools.
4. Fasten PVC firmly to socket with wire or rivets (or both). It is best to attach it temporarily with wire and to have the child try it before you fix it permanently.
5. Cover with sawdust and glue, or with resin-base casting bandage (very expensive) or with fiberglass and resin (also expensive).



Note: For a stronger, water-resistant limb, the socket can also be made with resin-based casting bandage. But this is also expensive.

Adjusting the PVC leg

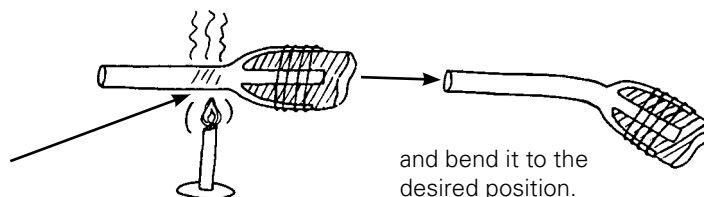After the child tries on the leg, if the angle seems wrong, you can gently heat the PVC near the stump,



and bend it to the desired position.

Lengthening the PVC leg to allow for growth of the child

1. Cut through PVC pipe.
2. Cut a strong hardwood rod or plastic tube that just fits inside the PVC.
3. Cut a ring of PVC pipe as long as the additional leg length needed. (Make it a little extra long to allow for more growth.)
4. Put the wood rod inside the PVC with the new ring of PVC.
5. Drill holes and screw or bolt leg together tightly.
6. When an even longer limb is needed, replace the short ring with a longer one.
7. If necessary, to make the joint stronger, cut a piece of PVC lengthwise, heat it, put it over the leg tube, and attach it tightly.

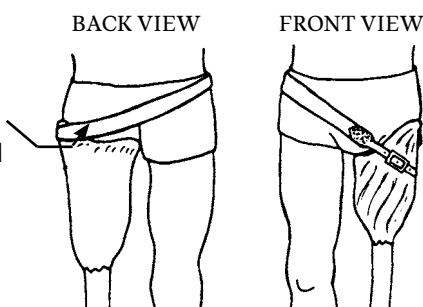
IMPORTANT: For both below-knee and above-knee limbs, try to line up the limb as well as possible so that its angle is similar to the other leg and feels right when the child stands and steps. Often this requires repeated tries and adjustments. Getting the limb to line up right is the key to successful limb fitting. It helps to learn this from someone skilled at fitting limbs.

ABOVE-KNEE ARTIFICIAL LIMBS

Children who are growing quickly need a low-cost limb that can be easily replaced or lengthened. Small children usually learn to walk well with a straight leg limb that does not have a knee joint.



1. A bamboo or PVC plastic tube above-knee limb can be made in much the same way as for the below-knee limb.

The top edge of the socket should be rounded to form a wide lip on the back, where the butt can sit. Weight bearing should be on the butt bone and over the entire stump—and not just on the end of the stump.




Note: In some countries, thin plastic cuffs the right shape for socket tops can be purchased in different sizes from orthopedic suppliers. They can be placed around the leg before casting and can be re-used. Ask for "prefabricated ischial weight-bearing cuffs."

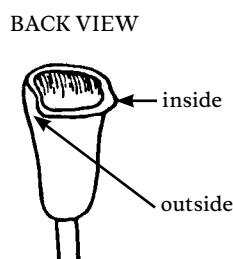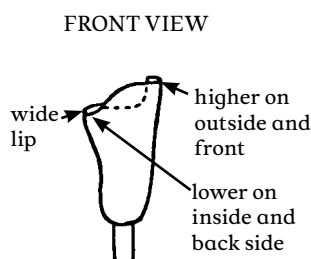Best shape for top of socket:




TOP VIEW

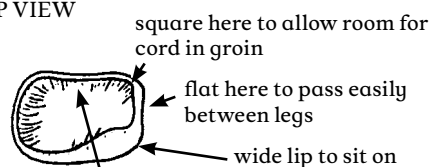


Make the top of the socket somewhat square. You can reshape the wet plaster a little while it is still on the stump, or just after you take it off.

2. A leather and metal rod limb (adapted from *Simple Prosthesis Manufacture*).

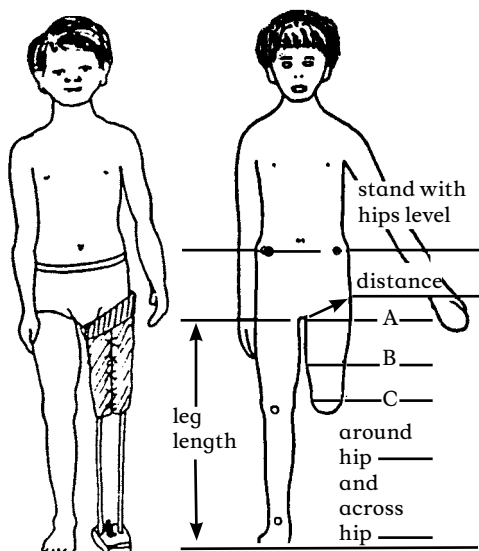MEASUREMENTS



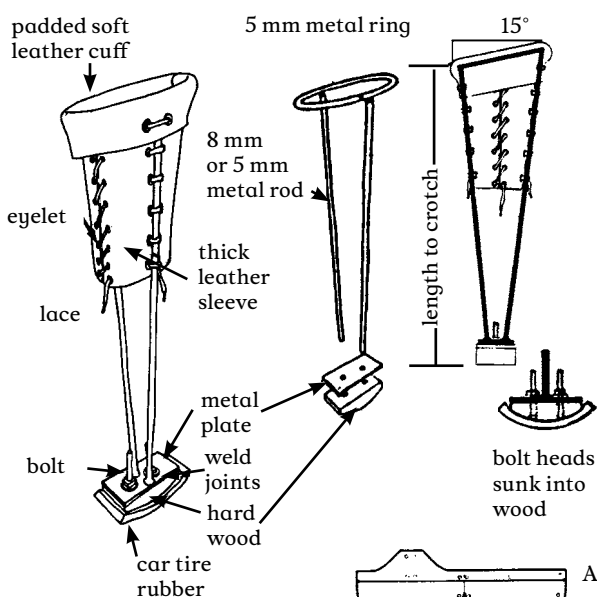NOTE: The socket is open at the bottom of the stump. This makes it cooler than the one above, and also allows for growth.




To make sleeve, draw on paper a design like this based on A, B, and C stump measurements. Then copy onto leather and cut it out.


Above-knee limb with knee joint (for older children and adults)

Adapted from the OHI Manual: *Simple Above-knee Prosthesis Manufacture*. (Manual out of print)

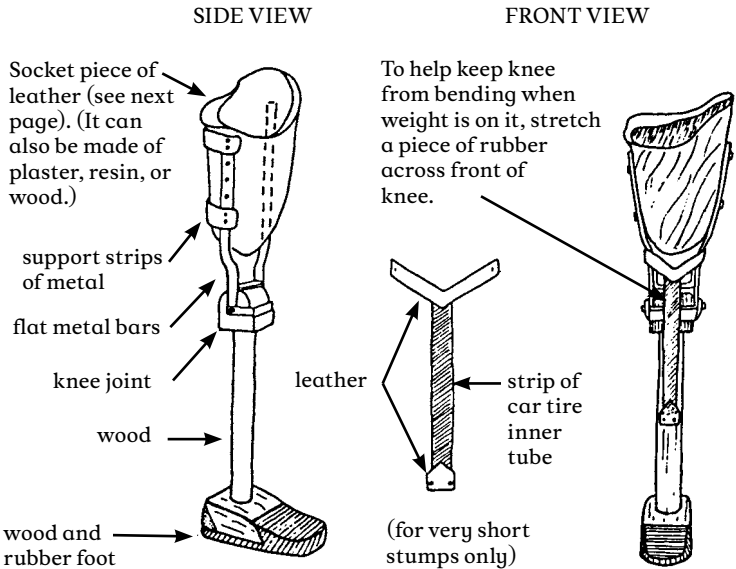NOTE: This is a simple design for a knee joint. Knee joints that work well are hard to make and you may need to experiment a lot. Perhaps you can re-use a joint from an old limb that is not being used.





Making the knee joint:

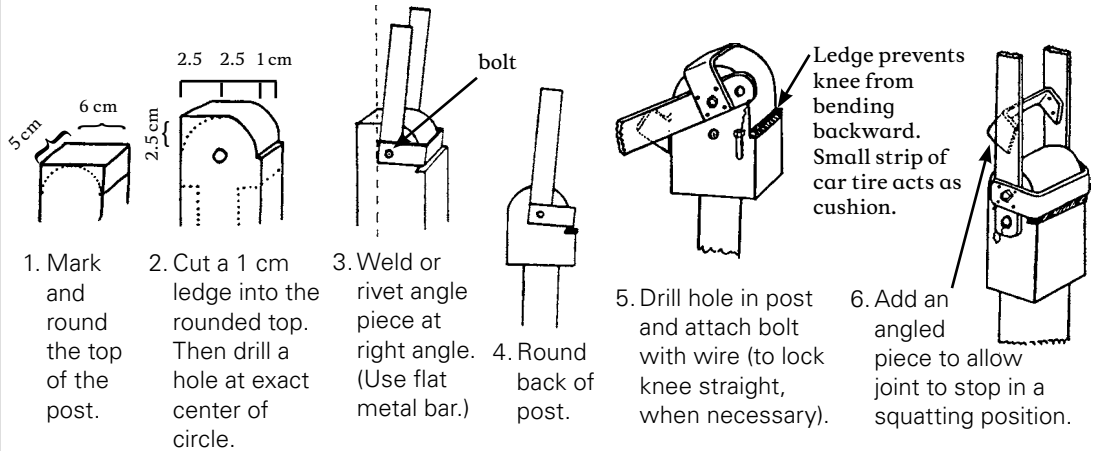

FEET

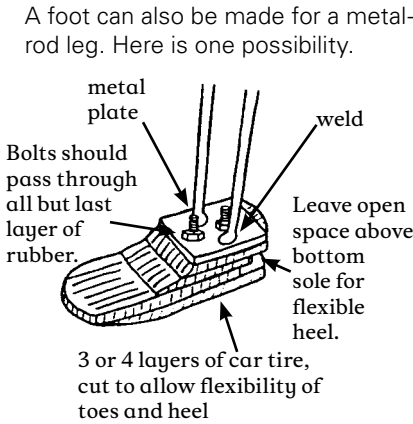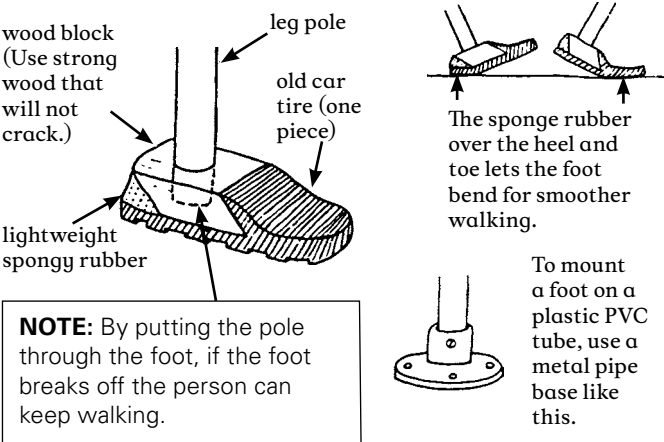
Putting feet on artificial legs makes them look better (with shoes, sandals, or boots). Also, the wide base helps prevent the leg from sinking into mud or sand. A well-made, flexible foot can make walking easier. Here are 2 possibilities.



OTHER WAYS OF MAKING ARTIFICIAL LIMBS

Wooden legs

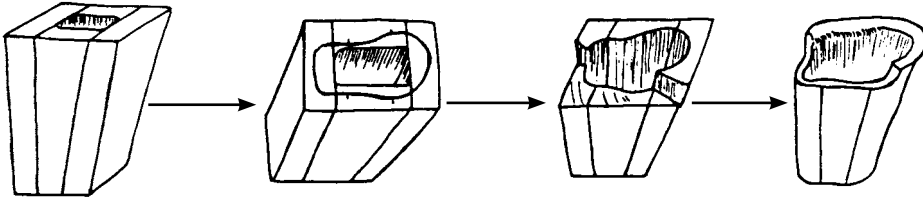The oldest, traditional way of making artificial limbs is to make the socket out of wood.




This is best learned from a skilled craftsman. Unfortunately, this is a skill that is difficult to learn from a book. A book that describes the method step-by-step is *Manual of Above-Knee Wood Socket Prosthetics* by Miles Anderson, John Bray, and Charles Hennessey. It has gone out of print, but you may still be able to find it. Unfortunately, the methods described are complex and require a lot of special equipment. However, perhaps they could be simplified. (We have not tried this method.)

Leather socket: Self-adjusting prosthesis

Adapted from the OHI Manual: *Simple Below-knee Prosthesis Manufacture*. (Manual out of print)


This method uses flat metal bars, a wood post, and a thick, firm leather socket. To form the socket, wet leather is stretched over a plaster mold of the stump.

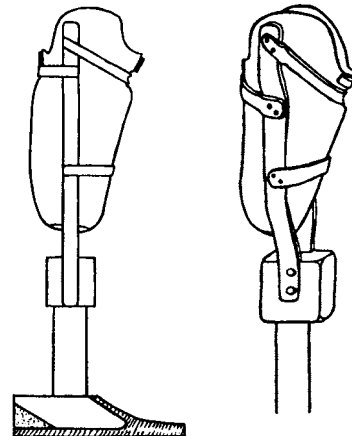A leather socket has several advantages. Leather is available almost everywhere, is more comfortable in hot weather, and can easily be adjusted to the stump as it becomes smaller. Also, leather is soft and easily takes the shape of the stump, and therefore self-corrects molding mistakes.



Stump protectors from old tires

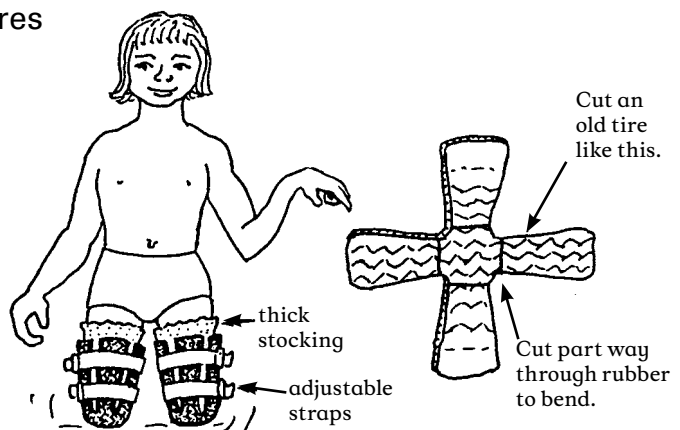
For a child with both legs amputated above the knee, short artificial limbs or even simple "stump protectors" may allow her to move about easier than long leg limbs.



Jaipur limb

The Jaipur limb was developed in Jaipur, India to meet the need for a limb that would (1) allow working barefoot in rice paddies, (2) look like a real bare foot, (3) bend at the foot in all directions enough so the person can squat easily and walk firmly on uneven ground, and (4) be low cost and quick to make.

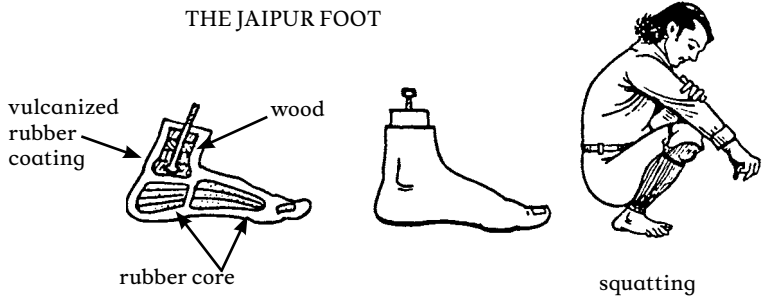The foot is made of wood and sponge rubber and then vulcanized (heat molded) with rubber, using a metal mold. The rubber gives the foot its life-like form and color and makes it strong and waterproof.



The limb is made of SHEET ALUMINUM.

1. Take measurements of full leg.

(Total cost for the below-knee limb can be less than US \$20.)

2. Mark a sheet of aluminum and cut out leg form.

3. Bend into tube.

4. Weld the seam.

5. Hammer tube into shape of leg. With repeated tests, fit it to the knee just right.

6. Attach knee strap and foot.

To make the Jaipur limb requires a lot of skill as well as special equipment. But once a shop is set up and persons trained, the limb can be made at very low cost, and fitted very quickly (one hour from the first measurements until the person walks away on his new limb). For instructions, contact Rehab Centre, SMS Medical College, Jaipur 302004, India.

Ideas for a limb-making shop. On p. 521 there is a description of the OHI prosthetics shop in Thailand, where workers with amputations make the bamboo and above-knee adjustable limbs shown in this chapter.

JAIPUR ABOVE-KNEE DESIGN

